BPL_TEST2_Perfusion - demo

In [1]: run -i BPL_TEST2_Perfusion_explore.py

This notebook explore perfusion cultivation in comparison with ordinary continuous cultivation (chemostat) and use comparable settings to earlier notebook. Further you see here examples of interaction with the simplified commands par(), init(), simu() etc as well as direct interaction with the FMU which is called "model" here. The last simulation is always available in the workspace and called "sim_res". Note that describe() brings mainly up from descriptive information from the Modelica code from the FMU but is complemented by some information given in the Python setup file.


```
Windows - run FMU pre-compiled JModelica 2.14
       Model for bioreactor has been setup. Key commands:
        - par()- change of parameters and initial values
        - init()

    change initial values only


       - simu()
                     - simulate and plot
       - newplot() - make a new plot
                   - show plot from previous simulation
        - show()
                    - display parameters and initial values from the last simulation
        - disp()
        - describe() - describe culture, broth, parameters, variables with values/units
       Note that both disp() and describe() takes values from the last simulation
       and the command process_diagram() brings up the main configuration
       Brief information about a command by help(), eg help(simu)
       Key system information is listed with the command system_info()
In [2]: %matplotlib inline
        plt.rcParams['figure.figsize'] = [25/2.54, 20/2.54]
```

No processDiagram.png file in the FMU, but try the file on disk.

In [3]: process_diagram()

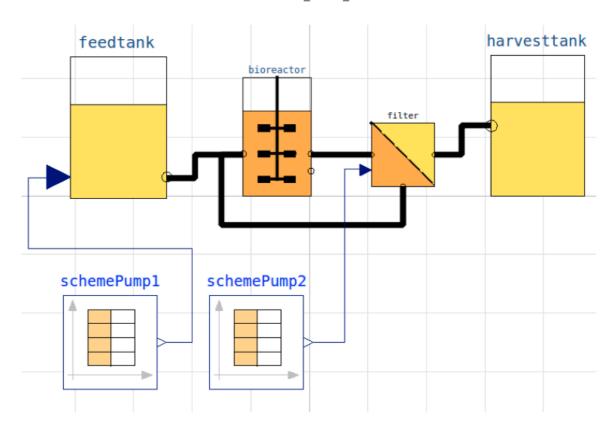

```
In [4]: # Process parameters used throughout
    par(Y=0.5, qSmax=0.75, Ks=0.1)  # Culture
    par(filter_eps=0.10, filter_alpha_X=0.02, filter_alpha_S=0.10)  # Filter
    par(S_in=30.0)  # Inlet subs
    init(V_start=1.0, VX_start=1.0)  # Process in
    eps = parDict['filter_eps']  # Pump schea
```

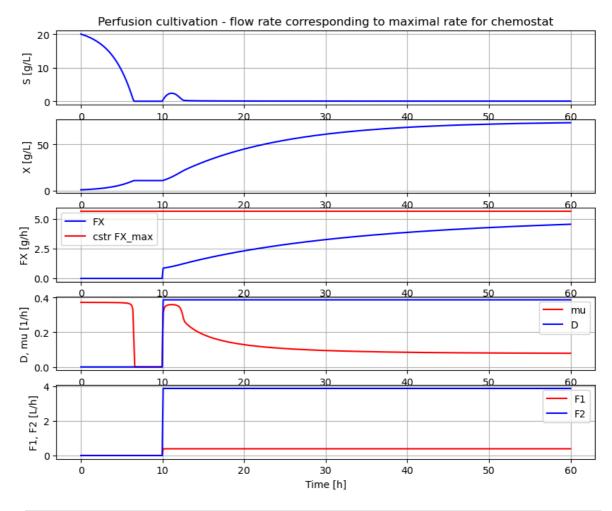


```
In [6]: # Concentration factor of the filter
    c=model.get('filter.retentate.c[1]')[0]/model.get('filter.inlet.c[1]')[0]
    print('Conc factor of perfusion filter =', np.round(c,3))
```

Conc factor of perfusion filter = 1.179

```
In [7]: c_data=sim_res['filter.retentate.c[1]']/sim_res['filter.inlet.c[1]']
print('Conc factor variation', np.round(min(c_data[151:]), 3),'to', np.round(max
```

Conc factor variation 1.179 to 1.179

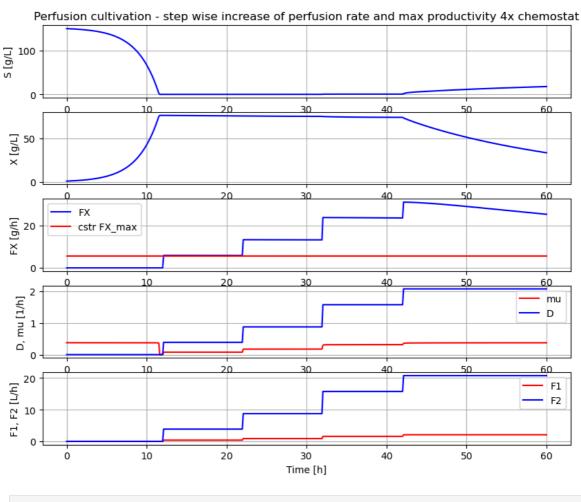

Conc factor variation 1.162 to 1.179

```
In [10]: describe('cstrProdMax')
```

Calculate from the model maximal chemostat productivity FX_max : 5.625 [g/h]

Out[11]: 23.6

Out[12]: 4.2

In [13]: # Finally we check the filter flow rates at time 40 hour - note the negative sig
model.get('filter.inlet.F')[0]

Out[13]: 15.74999999999998

In [14]: model.get('filter.filtrate.F')[0]

Out[14]: -1.575

```
In [15]: model.get('filter.retentate.F')[0]
Out[15]: -14.17499999999999
```

Summary

- The perfusion filter had a concentration factor of cells around 1.08 and re-cycling flow was set to a factor 10 higher than the perfusion rate and changed when perfusion rate was change to keep the ratio factor 10.
- The first simulation showed that by cell retention using perfusion filter the process could be run at a perfusion flow rate at the maximal flow rate possible for corresponding chemostat culture and cell concetration increased steadily.
- The second simulation showed that with a proper startup cell concentration, the cell
 concentration remained constant when perfusion rate increased in a similar way as
 what we see in a chemostat.
- The second simulation also showed that biomass productivity in this case was increased by a factor 4.2 compared to chemostat.
- If the perfusion rate increased to higher levels washout started but the decrase of cell concentration was slow.

Some of you who read this may have your perfusion experience with CHO-cultures. For such cultures the cell concentration do increase with increase of perfusion rate and there are understood reasons for that. But for this simplified process as well as microbial processes they typically keep cell concentration constant when flow rate is chaged, and that under quite wide conditions. I will try come back to this phenomena in a later notebook.

Appendix

 ${\bf System\ information}$

-OS: Windows
-Python: 3.10.14

-Scipy: not installed in the notebook

-PyFMI: 2.13.0

-FMU by: JModelica.org

-FMI: 2.0

-Type: FMUModelCS2

-Name: BPL_TEST2.Perfusion
-Generated: 2024-05-12T21:00:01

-MSL: 3.2.2 build 3

-Description: Bioprocess Library version 2.2.0

-Interaction: FMU-explore version 1.0.0

In []: