## BPL\_TEST2\_Perfusion script with FMPy

The key library FMPy is installed.

After the installation a small application BPL\_TEST2\_Perfusion is loaded and run. You can continue with this example if you like.

!lsb\_release -a # Actual VM Ubuntu version used by Google No LSB modules are available. Distributor ID: Ubuntu Description: Ubuntu 22.04.3 LTS Release: 22.04 Codename: jammy %env PYTHONPATH= → env: PYTH0NPATH= !wget https://repo.anaconda.com/miniconda/Miniconda3-py312\_24.3.0-0-Linux-x86\_64. !chmod +x Miniconda3-py312\_24.3.0-0-Linux-x86\_64.sh !bash ./Miniconda3-py312 24.3.0-0-Linux-x86 64.sh -b -f -p /usr/local import sys sys.path.append('/usr/local/lib/python3.12/site-packages/') → --2024-08-13 08:57:20-- https://repo.anaconda.com/miniconda/Miniconda3-py312 Resolving repo.anaconda.com (repo.anaconda.com)... 104.16.191.158, 104.16.32. Connecting to repo.anaconda.com (repo.anaconda.com)|104.16.191.158|:443... con HTTP request sent, awaiting response... 200 OK Length: 143351488 (137M) [application/octet-stream] Saving to: 'Miniconda3-py312\_24.3.0-0-Linux-x86\_64.sh' Miniconda3-py312\_24 100%[===========] 136.71M in 1.0s 131MB/s 2024-08-13 08:57:21 (131 MB/s) - 'Miniconda3-py312\_24.3.0-0-Linux-x86\_64.sh' : PREFIX=/usr/local Unpacking payload ... Installing base environment... Preparing transaction: ...working... done Executing transaction: ...working... done installation finished. !conda update -n base -c defaults conda --yes



BPL\_TEST2\_Perfusion\_fmpy\_colab.ipynb - Colab Upenss (-3.0.14 DIZ IID ע% שין [ששושש<י, יבור/5] 1.2 MB conda-24.7.1 0% 0/1 [00:00<?, ?it/s] certifi-2024.7.4 | 159 KB 0% 0/1 [00:00<?, ?it/s] | : ca-certificates-2024 | 127 KB 0% 0/1 [00:00<?, ?it/s] 1 : openssl-3.0.14 | 5.2 MB 0% 0.003006342237126712/1 [00:00<02:13 conda-24.7.1 1% 0.013060714305643354/1 [00:00<00:30 1 1.2 MB ca-certificates-2024 | 127 KB 13% 0.12647440251960723/1 [00:00<00:02, | : certifi-2024.7.4 I 159 KB 10% 0.10045740493212503/1 [00:00<00:03, 1 : openssl-3.0.14 | 5.2 MB 1: 47% 0.46598304675464036/1 [00:00<00:00, ca-certificates-2024 | 127 KB | : 100% 1.0/1 [00:00<00:00, 2.40it/s]

ca-certificates-2024 | 127 KB | : 100% 1.0/1 [00:00<00:00, 2.40it/s]

frozendict-2.4.2 | 36 KB | : 100% 1.0/1 [00:00<00:00, 1.97it/s]

frozendict-2.4.2 | 36 KB | : 100% 1.0/1 [00:00<00:00, 1.97it/s] certifi-2024.7.4 | 159 KB | : 100% 1.0/1 [00:00<00:00, 2.02it/s] certifi-2024.7.4 | 159 KB | : 100% 1.0/1 [00:00<00:00, 2.02it/s] conda-24.7.1 | 1.2 MB | : 100% 1.0/1 [00:01<00:00, 1.07s/it]

Preparing transaction: done Verifying transaction: done Executing transaction: done

!conda --version
!python --version

conda 24.7.1 Python 3.12.2 !conda install -c conda-forge fmpy --yes # Install the key package

 $\overline{\Rightarrow}$ 

executing transaction: done

```
#!conda install -c conda-forge matplotlib --yes
#!conda install -c conda-forge scipy --yes
#!conda install -c conda-forge openpyxl --yes
#!conda install -c conda-forge xlrd --yes
```

## Notes of BPL\_TEST2\_Perfusion

This notebook explore perfusion cultivation in comparison with ordinary continuous cultivation (chemostat) and use comparable settings to earlier notebook. Further you see here examples of interaction with the simplified commands par(), init(), simu() etc as well as direct interaction with the FMU which is called "model" here. The last simulation is always available in the workspace and called "sim\_res". Note that describe() brings mainly up from descriptive information from the Modelica code from the FMU but is complemented by some information given in the Python setup file.

Now specific installation run a simulation and notebook for that Start with connecting to Github. Then upload the two files:

- FMU BPL\_TEST2\_Perfusion\_linux\_om\_me.fmu
- Setup-file BPL\_TEST2\_Perfusion\_fmpy\_explore.py

```
%%bash
git clone https://github.com/janpeter19/BPL_TEST2_Perfusion
Cloning into 'BPL_TEST2_Perfusion'...
%cd BPL_TEST2_Perfusion
/content/BPL_TEST2_Perfusion
run -i BPL_TEST2_Perfusion_fmpy_explore.py
→ Linux – run FMU pre-comiled OpenModelica 1.23.0-dev
    Model for bioreactor has been setup. Key commands:
     - par()

    change of parameters and initial values

     - init()

    change initial values only

                   simulate and plot
     - simu()
     - newplot()
                   - make a new plot
     - show()
                   show plot from previous simulation
     - disp()

    display parameters and initial values from the last simulation
```

- describe() - describe culture, broth, parameters, variables with values/ur


Note that both disp() and describe() takes values from the last simulation and the command process\_diagram() brings up the main configuration

Brief information about a command by help(), eg help(simu)
Key system information is listed with the command system\_info()

%matplotlib inline
plt.rcParams['figure.figsize'] = [25/2.54, 20/2.54]

process\_diagram()


No processDiagram.png file in the FMU, but try the file on disk.

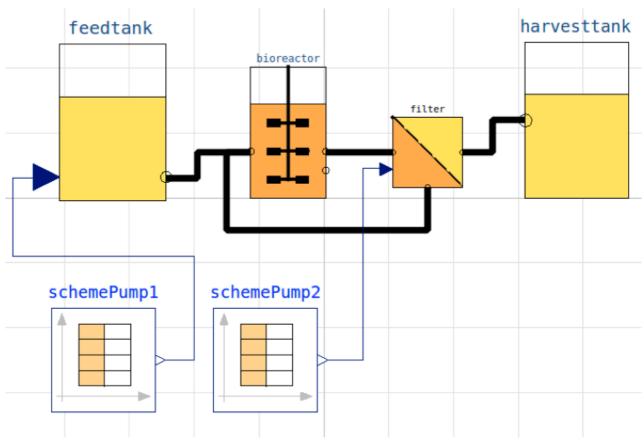

```
# Process parameters used throughout
par(Y=0.5, qSmax=0.75, Ks=0.1)  # Culture
par(filter_eps=0.10, filter_alpha_X=0.02, filter_alpha_S=0.10)  # Filter
par(S_in=30.0)  # Inlet subst
init(V_start=1.0, VX_start=1.0)  # Process ini
eps = parDict['filter_eps']  # Pump schedu
```

# Simulation of process with flow rate clot to wash-out for chemostat


newplot(title='Perfusion cultivation - flow rate corresponding to maximal rate fo simu(60)

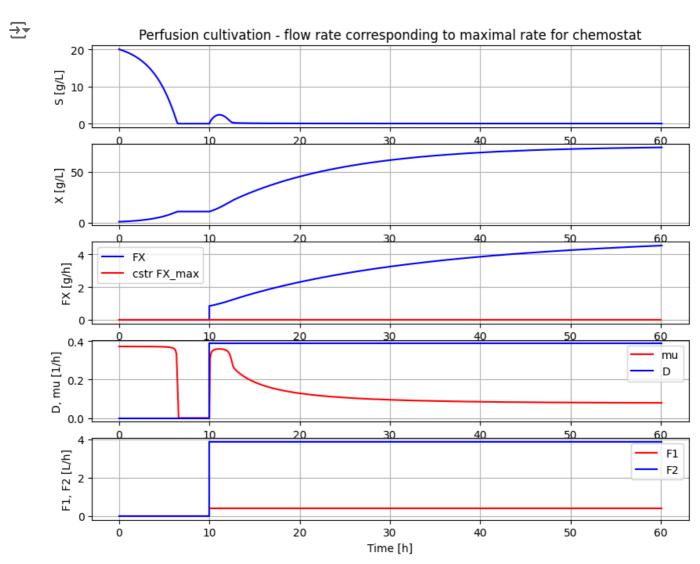


# Simulation of process with flow rate close to wash-out for chemostat

```
init(VS start=20)
                                                         # Process initial
par(pump1_t1=10, pump2_t1=10)
                                                         # Pump schedule - recycl
par(pump1_F1=2.5*0.155, pump2_F1=2.5*0.155/eps)
par(pump1_t2=940, pump2_t2=940, pump1_t3=950, pump2_t3=950, pump1_t4=960, pump2_t
newplot(title='Perfusion cultivation - flow rate corresponding to maximal rate fo
```

simu(30)simu(30,'cont')

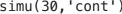

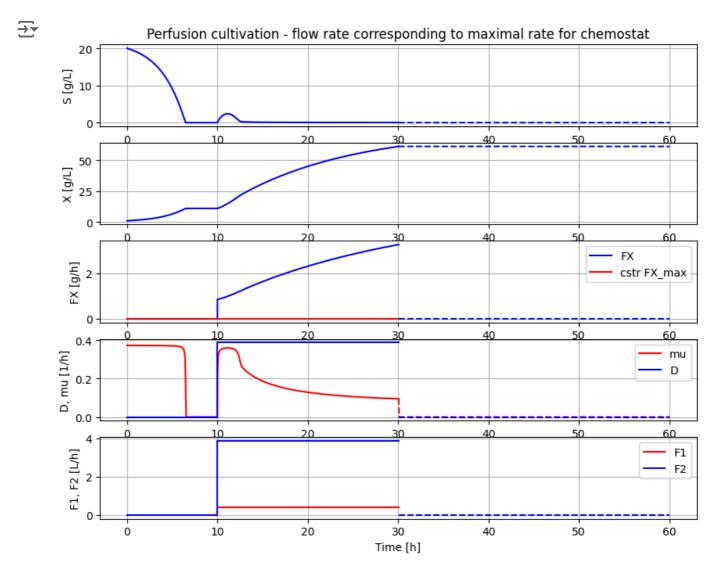


```
# Concentration factor of the filter
c=model_get('filter.retentate.c[1]')/model_get('filter.inlet.c[1]')
print('Conc factor of perfusion filter =', np.round(c,3))
```

Conc factor of perfusion filter = 1.186

```
c_data=sim_res['filter.retentate.c[1]']/sim_res['filter.inlet.c[1]']
print('Conc factor variation', np.round(min(c_data[151:]), 3), np.round(max(c_dat
Section Conc factor variation 1.186 1.186
# Simulation of process with step-wise increase of pefusion rate until wash-out.
# This means that re-circulation rate change at the same time as the perfusion ra
                                                          # Process initial varied
init(VS start=150)
par(pump1_t1=12, pump2_t1=12)
                                                          # Pump schedule - recycl
par(pump1 F1=2.5*0.155, pump2 F1=2.5*0.155/eps)
par(pump1_t2=22, pump2_t2=22)
par(pump1_F2=2.5*0.35, pump2_F2=2.5*0.35/eps)
par(pump1 t3=32, pump2 t3=32)
par(pump1_F3=2.5*0.63, pump2_F3=2.5*0.63/eps)
par(pump1_t4=42, pump2_t4=42)
par(pump1_F4=2.5*0.83, pump2_F4=2.5*0.83/eps)
```

newplot(title='Perfusion cultivation - step wise increase of perfusion rate and m simu(60)



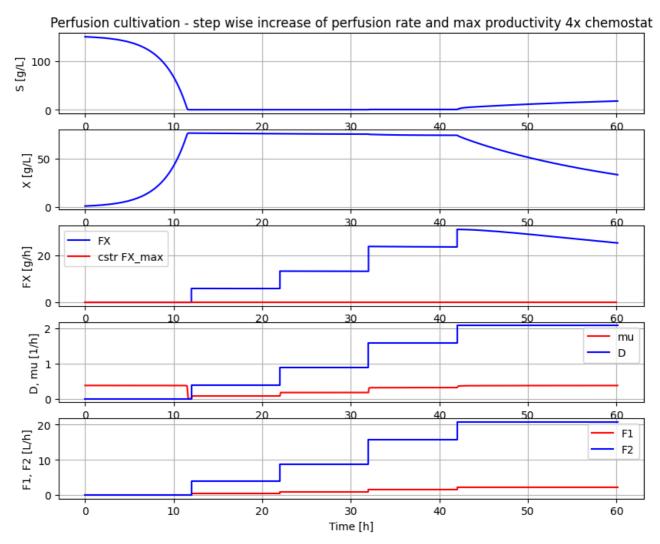

```
# Simulation without a plot and just to check typical values at high production r
#simu(40)
#c_data=sim_res['filter.retentate.c[1]']/sim_res['filter.inlet.c[1]']
#print('Conc factor variation', np.round(min(c_data[190:]), 3), 'to', np.round(ma
#describe('cstrProdMax')
```

# The maximal biomass productivity before washout is obtained aroudn 40 hours
np.round(model\_get('harvesttank.inlet.F')\*model\_get('harvesttank.inlet.c[1]'),1)

→ 25.2

# Thus perfusion (with this filter) brings a productivity improvement of about np.round(23.5/5.6,1)

```
→ 4.2
```

```
# Finally we check the filter flow rates at time 40 hour - note the negative sign model_get('filter.inlet.F')

→ 20.74999999999996

model_get('filter.filtrate.F')

→ -2.07499999999997

model_get('filter.retentate.F')

→ -18.674999999999997
```

## Summary

- The perfusion filter had a concentration factor of cells around 1.08 and re-cycling flow was set to a factor 10 higher than the perfusion rate and changed when perfusion rate was change to keep the ratio factor 10.
- The first simulation showed that by cell retention using perfusion filter the process could be run at a perfusion flow rate at the maximal flow rate possible for corresponding chemostat culture and cell concetration increased steadily.
- The second simulation showed that with a proper startup cell concentration, the cell
  concentration remained constant when perfusion rate increased in a similar way as what
  we see in a chemostat.
- The second simulation also showed that biomass productivity in this case was increased by a factor 4.2 compared to chemostat.
- If the perfusion rate increased to higher levels washout started but the decrase of cell concentration was slow.

Some of you who read this may have your perfusion experience with CHO-cultures. For such cultures the cell concentration do increase with increase of perfusion rate and there are understood reasons for that. But for this simplified process as well as microbial processes they typically keep cell concentration constant when flow rate is chaged, and that under quite wide conditions. I will try come back to this phenomena in a later notebook.

```
# List of components in the process setup and also a couple of other things like l: describe('parts')

['bioreactor', 'bioreactor.culture', 'D', 'feedtank', 'filter', 'harvesttank' describe('MSL')
```