BPL_TEST2_Perfusion script with FMPy

The key library FMPy is installed.

After the installation a small application BPL_TEST2_Perfusion is loaded and run. You can continue with this example if you like.

```
In [1]: !lsb_release -a # Actual VM Ubuntu version used by Google
       No LSB modules are available.
       Distributor ID: Ubuntu
                       Ubuntu 22.04.4 LTS
       Description:
       Release:
                       22.04
       Codename:
                       jammy
In [2]: !python --version
       Python 3.11.11
In [3]: !pip install fmpy
       Collecting fmpy
         Downloading FMPy-0.3.22-py3-none-any.whl.metadata (1.9 kB)
       Requirement already satisfied: attrs in /usr/local/lib/python3.11/dist-packages (fro
       m fmpy) (25.1.0)
       Requirement already satisfied: Jinja2 in /usr/local/lib/python3.11/dist-packages (fr
       om fmpy) (3.1.5)
       Collecting lark (from fmpy)
         Downloading lark-1.2.2-py3-none-any.whl.metadata (1.8 kB)
       Requirement already satisfied: lxml in /usr/local/lib/python3.11/dist-packages (from
       fmpy) (5.3.1)
       Requirement already satisfied: msgpack in /usr/local/lib/python3.11/dist-packages (f
       rom fmpy) (1.1.0)
       Requirement already satisfied: numpy in /usr/local/lib/python3.11/dist-packages (fro
       m fmpy) (1.26.4)
       Requirement already satisfied: MarkupSafe>=2.0 in /usr/local/lib/python3.11/dist-pac
       kages (from Jinja2->fmpy) (3.0.2)
       Downloading FMPy-0.3.22-py3-none-any.whl (4.9 MB)
                                                 - 4.9/4.9 MB 21.2 MB/s eta 0:00:00
       Downloading lark-1.2.2-py3-none-any.whl (111 kB)
                                                 - 111.0/111.0 kB 4.9 MB/s eta 0:00:00
       Installing collected packages: lark, fmpy
       Successfully installed fmpy-0.3.22 lark-1.2.2
```

Notes of BPL_TEST2_Perfusion

This notebook explore perfusion cultivation in comparison with ordinary continuous cultivation (chemostat) and use comparable settings to earlier notebook. Further you see here examples of interaction with the simplified commands par(), init(), simu() etc as well as direct interaction with the FMU which is called "model" here. The last simulation is always


available in the workspace and called "sim_res". Note that describe() brings mainly up from descriptive information from the Modelica code from the FMU but is complemented by some information given in the Python setup file.

Now specific installation run a simulation and notebook for that Start with connecting to Github. Then upload the two files:

- FMU BPL_TEST2_Perfusion_linux_om_me.fmu
- Setup-file BPL_TEST2_Perfusion_fmpy_explore.py

```
In [4]: | %%bash
        git clone https://github.com/janpeter19/BPL_TEST2_Perfusion
      Cloning into 'BPL_TEST2_Perfusion'...
In [5]: %cd BPL_TEST2_Perfusion
       /content/BPL_TEST2_Perfusion
In [6]: run -i BPL_TEST2_Perfusion_fmpy_explore.py
       Linux - run FMU pre-comiled OpenModelica
       Model for the process has been setup. Key commands:
        par()change of parameters and initial values
       init()change initial values onlysimu()simulate and plot
        - newplot() - make a new plot
        show()show plot from previous simulation
                    - display parameters and initial values from the last simulation
        - disp()
        - describe() - describe culture, broth, parameters, variables with values/units
       Note that both disp() and describe() takes values from the last simulation
       and the command process_diagram() brings up the main configuration
       Brief information about a command by help(), eg help(simu)
       Key system information is listed with the command system_info()
In [7]: %matplotlib inline
        plt.rcParams['figure.figsize'] = [25/2.54, 20/2.54]
In [8]: process_diagram()
```


No processDiagram.png file in the FMU, but try the file on disk.

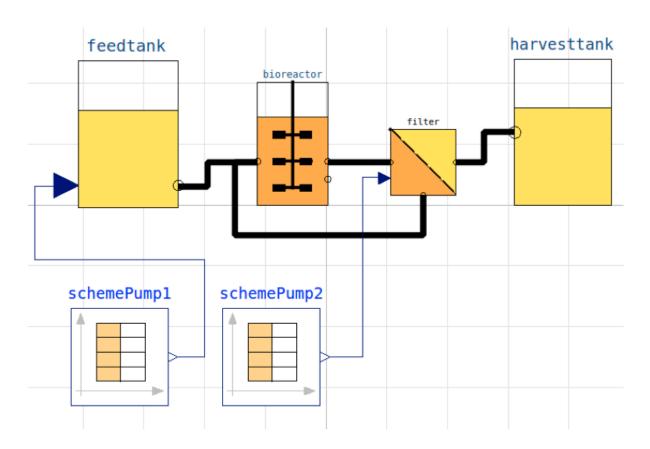

```
In [9]: # Process parameters used throughout
    par(Y=0.5, qSmax=0.75, Ks=0.1)  # Culture
    par(filter_eps=0.10, filter_alpha_X=0.02, filter_alpha_S=0.10)  # Filter
    par(S_in=30.0)  # Inlet substra
    init(V_start=1.0, VX_start=1.0)  # Process initi
    eps = parDict['filter_eps']  # Pump schedule
```

```
In [10]: # Simulation of process with flow rate clot to wash-out for chemostat


init(VS_start=20)  # Process initial
  par(pump1_t1=10, pump2_t1=10)  # Pump schedule - recycle
  par(pump1_F1=2.5*0.155, pump2_F1=2.5*0.155/eps)
  par(pump1_t2=940, pump2_t2=940, pump1_t3=950, pump2_t3=950, pump1_t4=960, pump2_t4=
  newplot(title='Perfusion cultivation - flow rate corresponding to maximal rate for simu(60)
```

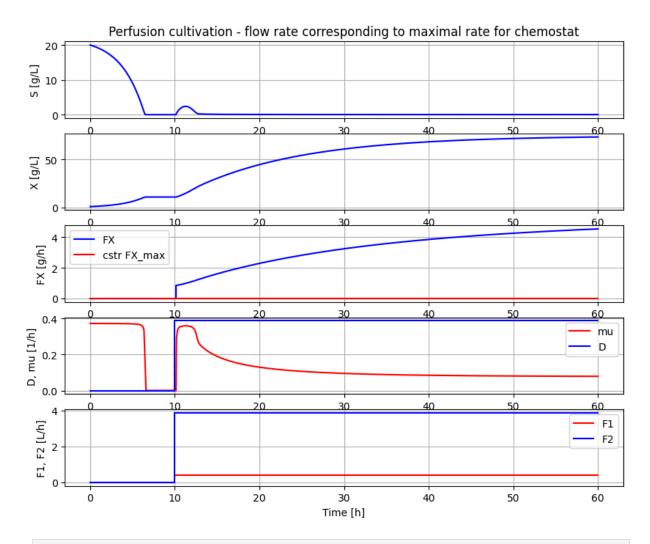

```
In [11]: # Simulation of process with flow rate close to wash-out for chemostat

init(VS_start=20)  # Process initial
    par(pump1_t1=10, pump2_t1=10)  # Pump schedule - recycle
    par(pump1_F1=2.5*0.155, pump2_F1=2.5*0.155/eps)
    par(pump1_t2=940, pump2_t2=940, pump1_t3=950, pump2_t3=950, pump1_t4=960, pump2_t4=

newplot(title='Perfusion cultivation - flow rate corresponding to maximal rate for
    simu(30)
    simu(30,'cont')
```

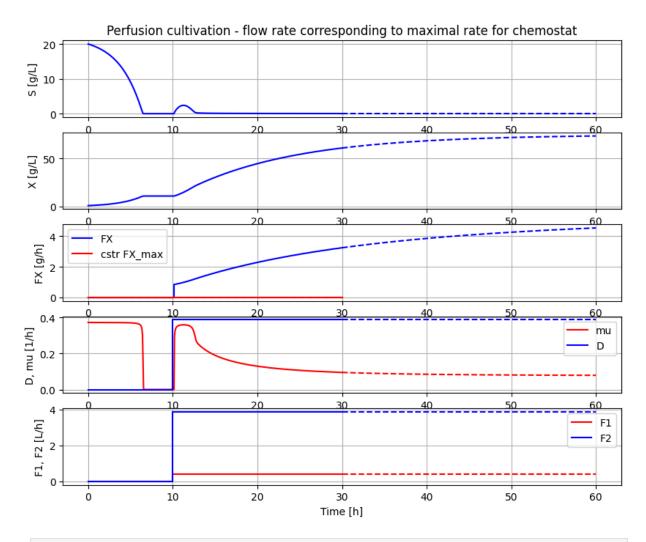


```
In [12]: # Concentration factor of the filter
    c=model_get('filter.retentate.c[1]')/model_get('filter.inlet.c[1]')
    print('Conc factor of perfusion filter =', np.round(c,3))
```

Conc factor of perfusion filter = 1.179

```
In [13]: c_data=sim_res['filter.retentate.c[1]']/sim_res['filter.inlet.c[1]']
    print('Conc factor variation', np.round(min(c_data[151:]), 3), np.round(max(c_data[
```

Conc factor variation 1.179 1.179

newplot(title='Perfusion cultivation - step wise increase of perfusion rate and max simu(60)

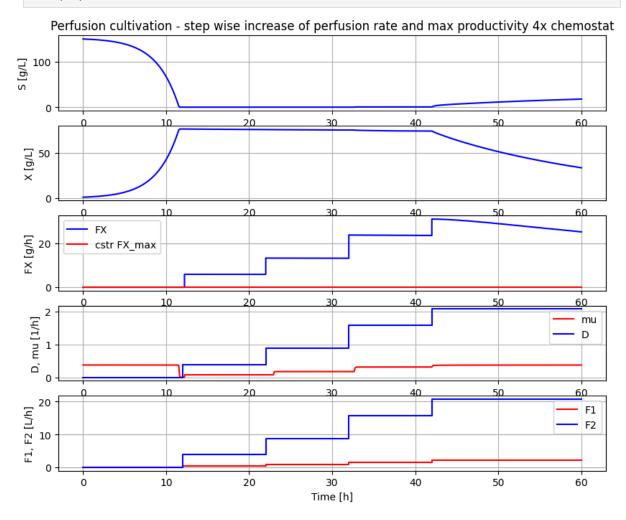

```
In [15]: # Simulation without a plot and just to check typical values at high production rat
#simu(40)
#c_data=sim_res['filter.retentate.c[1]']/sim_res['filter.inlet.c[1]']
#print('Conc factor variation', np.round(min(c_data[190:]), 3), 'to', np.round(max(
```

In [16]: #describe('cstrProdMax')

Out[17]: 25.2

Out[18]: 4.2

In [19]: # Finally we check the filter flow rates at time 40 hour - note the negative sign f
model_get('filter.inlet.F')

Out[19]: 20.74999999999996

```
In [20]: model_get('filter.filtrate.F')
Out[20]: -2.07499999999997
In [21]: model_get('filter.retentate.F')
Out[21]: -18.67499999999997
```

Summary

- The perfusion filter had a concentration factor of cells around 1.08 and re-cycling flow was set to a factor 10 higher than the perfusion rate and changed when perfusion rate was change to keep the ratio factor 10.
- The first simulation showed that by cell retention using perfusion filter the process could be run at a perfusion flow rate at the maximal flow rate possible for corresponding chemostat culture and cell concetration increased steadily.
- The second simulation showed that with a proper startup cell concentration, the cell concentration remained constant when perfusion rate increased in a similar way as what we see in a chemostat.
- The second simulation also showed that biomass productivity in this case was increased by a factor 4.2 compared to chemostat.
- If the perfusion rate increased to higher levels washout started but the decrase of cell concentration was slow.

Some of you who read this may have your perfusion experience with CHO-cultures. For such cultures the cell concentration do increase with increase of perfusion rate and there are understood reasons for that. But for this simplified process as well as microbial processes they typically keep cell concentration constant when flow rate is chaged, and that under quite wide conditions. I will try come back to this phenomena in a later notebook.

```
In [22]: # List of components in the process setup and also a couple of other things like li
    describe('parts')

['bioreactor', 'bioreactor.culture', 'D', 'feedtank', 'filter', 'harvesttank', 'sche
    mePump1', 'schemePump2']

In [23]: describe('MSL')

MSL: 3.2.3 - used components: RealInput, RealOutput, CombiTimeTable, Types

In [24]: system_info()
```

```
System information
```

-OS: Linux

-Python: 3.11.11

-Scipy: not installed in the notebook

-FMPy: 0.3.22

-FMU by: OpenModelica Compiler OpenModelica 1.25.0~dev-133-ga5470be

-FMI: 2.0 -Type: ME

-Name: BPL.Examples_TEST2.Perfusion -Generated: 2024-11-06T21:37:58Z

-MSL: 3.2.3

-Description: Bioprocess Library version 2.3.0 -Interaction: FMU-explore for FMPy version 1.0.1

In [24]: